Folding Massage Table
Table de Massage Pliante

HB88516

EN

Before You Start

- ⚠ Please read all instructions carefully.
- ⚠ Retain instructions for future reference.
- ⚠ Separate and count all parts and hardware.
- ⚠ Read through each step carefully and follow the proper order.
- ⚠ We recommend that, where possible, all items are assembled near to the area in which they will be placed in use, to avoid moving the product unnecessarily once assembled.
- ⚠ Always place the product on a flat, steady and stable surface.
- ⚠ Keep all small parts and packaging materials for this product away from babies and children as they potentially pose a serious choking hazard.

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
EN

Precautionary Instructions


1. Maintain the table periodically. Make sure all screws are secure and the cable is not deteriorated or frayed. Make sure there is no cracks on the wooden structure.
 2. Do not overload the table beyond the recommended weight limit. This could result in serious injury to the user. This product is tested to hold a max working load of 441 lbs (200 kg).
 3. Do not use alcohol based cleaners when cleaning the PVC. This could cause damage to the PVC material. Use warm soapy water and wipe dry.
 4. Make sure when using the table that all four legs are at the same level. Using the legs at different levels can cause the table to be unstable and cause injury.
 5. Do not allow children to assemble or disassemble the massage table for you unless they are under supervision.
 6. Be careful when using sharp objects as they can pierce the covering. Keep the table's covering away from being contacted by the sharp items.
- Warning:** The loading capacity for the headrest cradle: 10.0kgs.

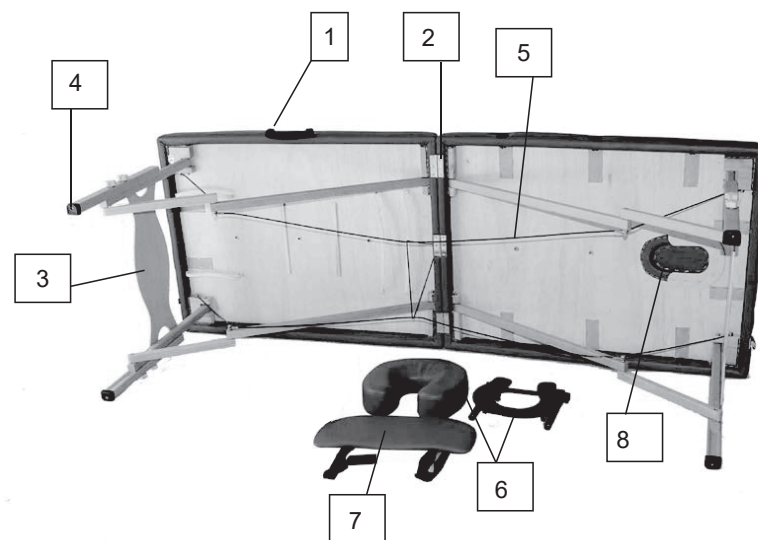- | | | |
|-----------------------|------------------|--------------|
| 1. Handle | 2. Copper button | 3. End plate |
| 4. Rubber foot | 5. Cable | 6. Headrest |
| 7. Hanged armrest pad | 8. Face hole | |

Setting Up Your Table

Steps are as below:

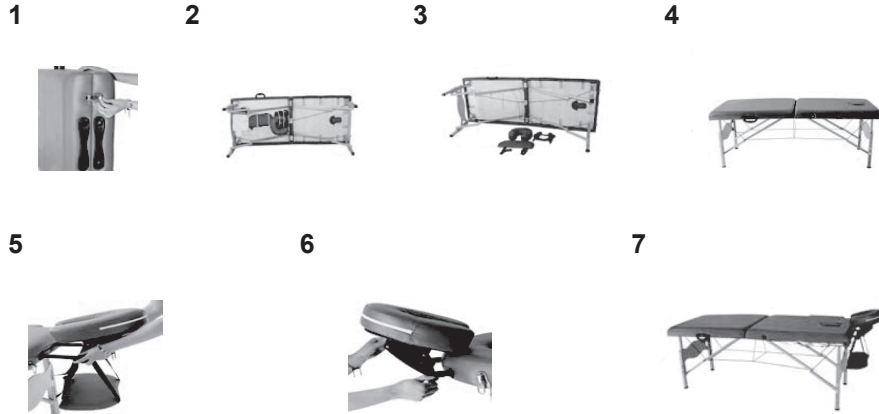1. Release the buckle.
2. Lay the table on its side with four rubber feet on the floor, open up the table with four legs extending as full as possible.
3. Take out all accessories.
4. Lift the table up until all legs are on the floor.
5. Put accessories into the corresponding outlets.
6. Release the side knob of the headrest cradle, and then press and fix it once the headrest is adjusted to the required position.
7. Now the table is ready for your use.

Re: The height of the table is adjustable. Turn off the knobs, take out the movable leg to adjust the height by choosing the holes on the movable leg, and put back onto the leg and then fasten the knobs firmly. Four legs should be done in the same way one by one and make sure all legs are adjusted into the same height level.

Disassemble and store the table

1. Lift up the backrest to the max and put down the backrest, lock the buckle.
2. Take out all accessories.
3. Lay the table flat with four rubber feet on the ground.
4. Put all accessories under the elastic strap under the top veneer.
5. Fold the legs on each end at angle inwards by pressing the corresponding support arms.
6. Fold the table up and lock the buckle.
7. Carry the table like suitcase by grasping the handle.

Warning: If the frame/legs/table structure is damaged in anyway, you should not use this table.

Using this table when damaged may cause harm or serious injury to you or your patient. You should seek repair immediately.

Warning card for the cables

While opening the massage table, please unfold the cables in a proper way, otherwise, the legs of the table will be broken and dangerous.

The right way to unfold the cables is as below:

The long and the short cables should be on a line.

The wrong way to unfold the cables is as below:

It is dangerous to wrap the cable around the wood tube.

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EN

Return / Damage Claim Instructions

- ⚠ **DO NOT discard the box / original packaging.**
In case a return is required, the item must be returned in original box. Without this your return will not be accepted.
- ⚠ **Take a photo of the box markings.**
A photo of the markings (text) on the side of the box is required in case a part is needed for replacement. This helps our staff identify your product number to ensure you receive the correct parts.
- ⚠ **Take a photo of the damaged part (if applicable).**
A photo of the damage is always required to file a claim and get your replacement or refund processed quickly. Please make sure you have the box even if it is damaged.
- ⚠ **Send us an email with the images requested.**
Email us directly from marketplace where your item was purchased with the attached images and a description of your claim.

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