

AP540022H



ALL ILLUSTRATIONS OF THE PRODUCT ARE FOR REFERENCE ONLY.
COMPONENTS SUPPLIED WITH INDIVIDUAL PRODUCTS MAY DIFFER.



ALWAYS READ THROUGH AND FULLY UNDERSTAND THE INSTRUCTION LEAFLET BEFORE STARTING THE INSTALLATION OF THIS PRODUCT.

ALWAYS CHECK THE PACKAGING AND ENSURE YOU HAVE ALL FIXINGS AND COMPONENT PARTS DETAILED IN THIS INSTRUCTION LEAFLET BEFORE RECYCLING OR DISPOSING OF THE PACKAGING.

ALWAYS ENSURE YOU HAVE ALL THE FIXINGS AND COMPONENT PARTS BEFORE ATTEMPTING TO ASSEMBLE OR FIT THIS PRODUCT.

Thank you for purchasing this Croydex product. Follow the fixing and care instructions in this manual to ensure that your product is fitted safely and that it retains its high quality finish.

RETAIN THIS MANUAL FOR FUTURE REFERENCE.

ADDRESS (UK)

Croydex
Central Way
Andover, Hampshire
SP10 5AW, UK

ADDRESS (EU)

MIL, Unit 8
Purcellsinch Ind Est
Dublin Road, Kilkenny
R95 HP71, Ireland


1 - INTRODUCTION

SAFETY & PRODUCT CARE


1. DO NOT use power tools to assemble this product.
2. DO NOT exceed the maximum user weight indicated in the fitting instructions.
3. ALWAYS check the legs are level and locked in place and the arm rests and back rests are secure.
4. Only use on surfaces free from obstacles.
6. ALWAYS assemble in accordance with the instructions provided.
7. ALWAYS make sure all four feet are in contact with the supporting surface at all times.
8. DO NOT modify the shower chair and its components in any way.
9. DO NOT use the shower chair if damaged or altered in any way.
10. REMOVE the shower chair after use.
12. DO NOT stand or kneel on the shower chair.
13. DO NOT use abrasive or chemical cleaners as these may damage the product. To retain the best quality finish, clean the product regularly with a soft damp cloth.

PRODUCT FITTING

IMPORTANT: Before attempting to install or use this product, check that it will fit and function in the desired location.

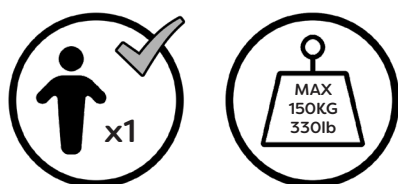


COMPONENTS (Not to scale)

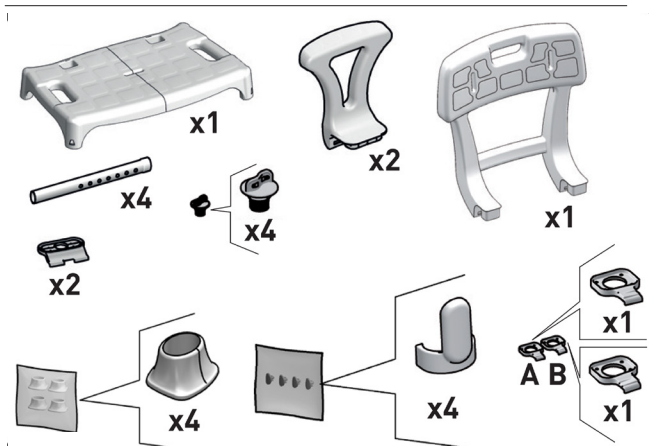

2 - INSTALLATION


STEP 1

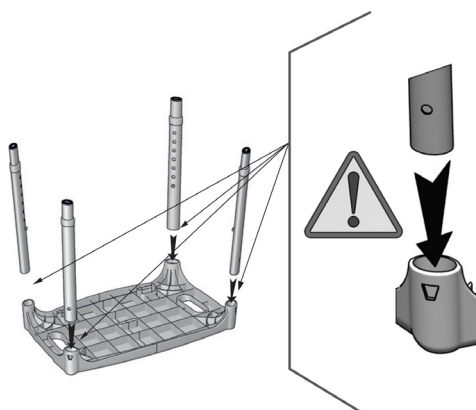

Firmly press all four adjustable legs into the bottom of the seat, ensuring that the adjustment holes are facing inwards and the fixing holes on the legs are aligned with the fixing holes in the seat.

Note: Ensure the thicker end of the leg is slotted into the seat to ensure the fixing holes line up and there is no movement.


STEP 2

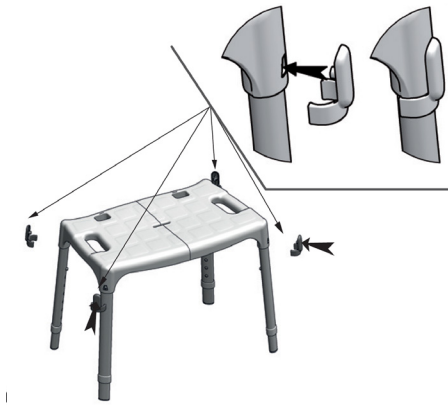

Press the seat plugs into the fixing holes to secure the adjustable legs to the seat.


Note: Ensure the seat plugs are aligned and orientated correctly.

STEP 3

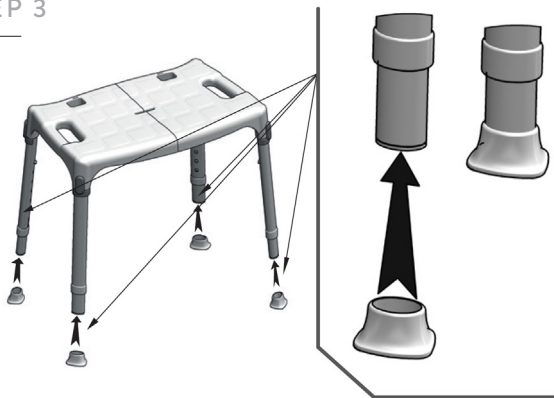

Firmly press the feet onto the bottom of each adjustable leg. Ensure the feet are angled outwards.


STEP 4

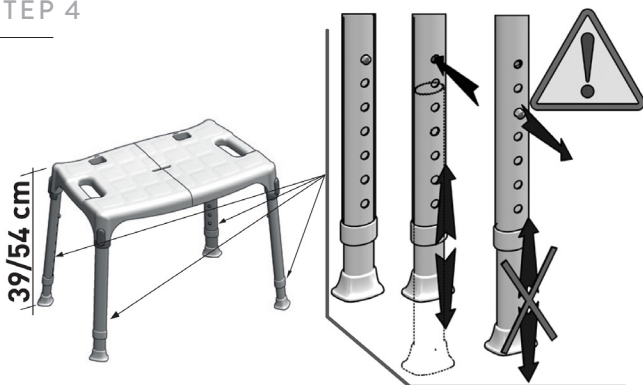

The height of the chair can be adjusted by depressing the locking button on each leg then sliding the thinner leg section up or down. Ensure that the locking button is relocated into one of the holes and check that the adjustable leg is locked in place before use.


YOU MUST ensure that all legs are fully secured and are at the same height before use.

STEP 5

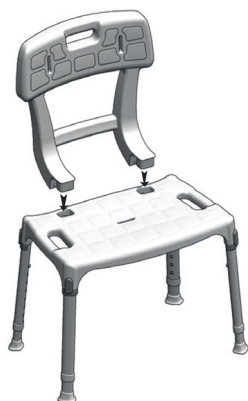

Slot the back rest into the two holes at the back of the chair.

STEP 6

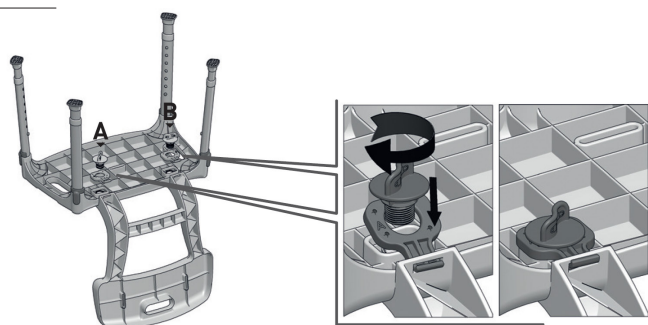


To secure the back rest in position, insert the two tabs on the bracket through the slots in the back rest. Push both brackets down fully and insert the back rest screws. Turn clockwise to secure the back rest in place. Ensure both screws are tightly fastened.

STEP 7



Slot both arm rests into the sides of the chair.

STEP 8

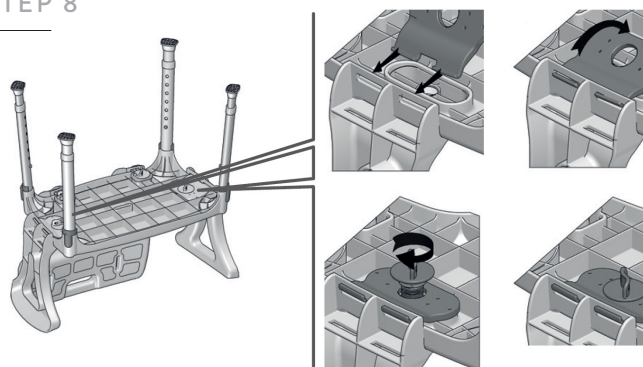

To secure the arm rest in position, insert the two tabs on the bracket through the slots in the arm rest. Push the bracket down fully and insert arm rest screw. Turn clockwise to secure the arm rest in place. Repeat this for the other arm rest. Ensure both screws are fastened tightly.

STEP 9



MAX USER
WEIGHT:
150 kg
330 lb

The chair is fully assembled and ready to be used.

Note: Ensure you have read the safety and product care section before use.